

# Respireo Primo F

simple & convenient facial mask



## We are all different

**The Respireo Primo series is a range of masks for ventilation that has been specifically designed to fit each need of the patients and adapt to each morphology.**

Simple & convenient, Respireo Primo F is a full face mask with adjustable dial and high ease of use. It comes in vented and non-vented and has been specifically designed to ensure your patients' maximum comfort.

It is also recommended for use with patients who mouth breath, have high pressures or are struggling to cope with therapy. The intuitive fastening clip and simple assembly make fitting this mask an easy task.

Respireo  
Primo **F**

# Respireo Primo F Mask




## Dual wall cushion

Ensuring stability and air-tightness with a great comfort thanks to the thin external lip.

Easy assembly with only 2 components.

## 'V' shaped clips


Flared clip system for an intuitive fitting and a quick release system if need be.

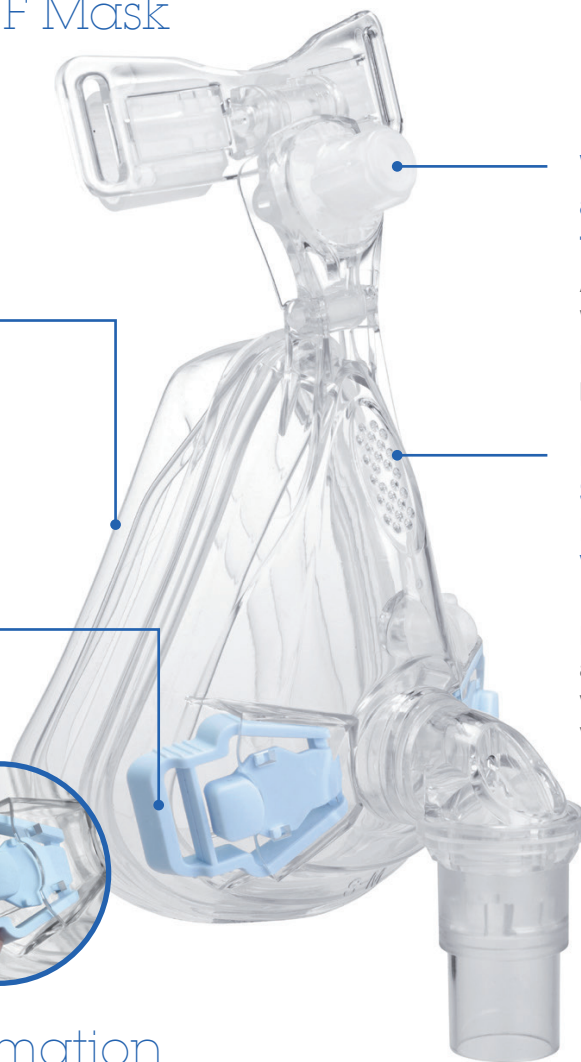

## Well distributed and adjustable forehead support

Adjusting dial to find the perfect fitting with the forehead inclination.

Large silicone pad to distribute the pressure points.

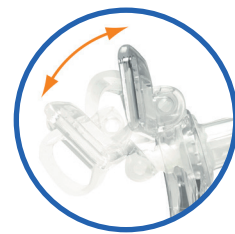

## Diffused and Silent Exhalation Valve

Respireo Primo F is also available in **non-vented version** for patients with valve ventilation.



## Ordering information

Product information	Size	Part number	Manufacturing reference
<b>Respireo Primo F</b>			
Respireo Primo F	Small	<b>KM217900</b>	(Manufacturer reference 472530*)
Respireo Primo F	Medium	<b>KM218000</b>	(Manufacturer reference 472531*)
Respireo Primo F	Large	<b>KM218100</b>	(Manufacturer reference 472532*)
<b>Respireo Primo F Non-Vented</b>			
Respireo Primo F Non-Vented	Small	<b>KM218200</b>	(Manufacturer reference 472533*)
Respireo Primo F Non-Vented	Medium	<b>KM218300</b>	(Manufacturer reference 472534*)
Respireo Primo F Non-Vented	Large	<b>KM218400</b>	(Manufacturer reference 472535*)
<b>Respireo Primo F Cushion</b>			
Respireo Primo F Cushion (pk 10)	Small	<b>KM218500</b>	(Manufacturer reference 472615)
Respireo Primo F Cushion (pk 10)	Medium	<b>KM218600</b>	(Manufacturer reference 472616)
Respireo Primo F Cushion (pk 10)	Large	<b>KM218700</b>	(Manufacturer reference 472617)
<b>Others</b>			
Headgear Respireo	Standard	<b>KM224000</b>	(Manufacturer reference 472456*)
Headgear Respireo	Large	<b>KM224100</b>	(Manufacturer reference 472457*)
Respireo Primo F Sizing Gauge (pk 20)		<b>KM220400</b>	(Manufacturer reference 472432)

(\*) Medical device class IIa 0051.

For more information please contact us at:

Air Liquide Healthcare, Stratus House, Unit 1 IDA College Park, Blanchardstown Road North, Blanchardstown, Dublin 15, D15 PEC4.

Tel: +353 (0)1 809 1800 Freephone: 1800 24 02 02 (ROI only) 0800 328 5875 (NI only)

Fax: +353 (0)1 829 3966 Email: [healthie@airliquide.ie](mailto:healthie@airliquide.ie) Website: [www.airliquidehealthcare.ie](http://www.airliquidehealthcare.ie)